


# Instructions Manual

## Beauty instrument

WKB017

DIAMOND DERMABRASION

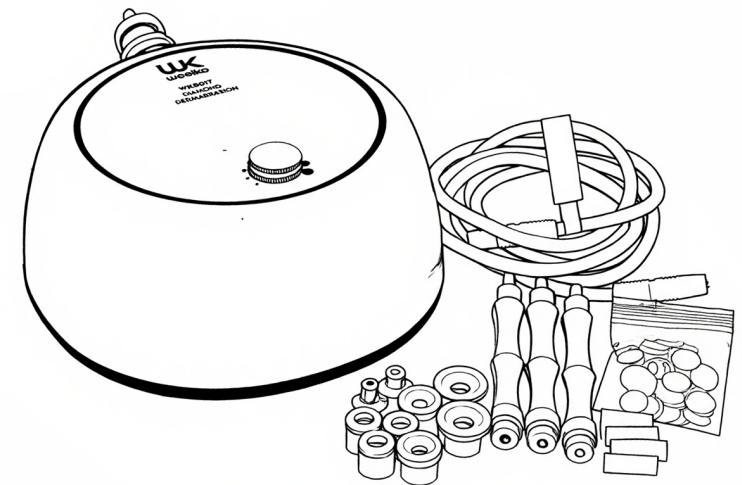

**WK**  
Imported by  
**WEEKO BARCELONA SL**

Pol. Ind. La Vall d'an  
Camí de Garreta, 92  
08600 Berga  
Barcelona - Spain  
T. +34 93 824 90 20  
E. info@weelko.com

VERSION 1.0

Thank you sincerely for choosing our product. This manual is compiled to help users operate this instrument effectively. We hope our company's products and services will bring you a distinguished experience.

This user manual contains all general information about the instrument. It serves as the starting point for new users to begin using the instrument, as it provides a general introduction to the system, explains how it works, and details the instrument's installation and various operation methods.

### Usage Precautions

1. Before turning on the machine, be sure to read the user manual carefully. Operators must be familiar with the user manual and operation methods, and the instrument must be managed by a designated person;
2. Do not use a power strip to avoid fire hazards;
3. When replacing the fuse, always turn off the power switch to prevent fire and electric shock;
4. Must use a power supply of 220V±22V, 50HZ±1HZ and an effectively grounded socket;
5. Using a damaged or modified power cord may cause fire and electric shock hazards;
6. If the instrument is loose, falling, or damaged, please contact the company in a timely manner;
7. Do not use the instrument in unstable environments such as tilting, vibrating, or impacting conditions;
8. No water or debris is allowed to enter the interior of the instrument. In case of foreign objects entering, please contact the company promptly;
9. When moving the instrument, the power plug must be pulled out first;
10. When moving the instrument, try to avoid vibration;
11. The instrument should be placed on a stable surface, with at least 5cm of space left around it to ensure air circulation and facilitate heat dissipation of the instrument;
12. Do not cover the instrument with anything to prevent blockage of the air vents;
13. Avoid using the instrument in the following environments: direct sunlight; high humidity environments; near water; vibrating or tilting places; strong magnetic field environments; places where electromagnetic waves or impulse voltages occur; places where chemicals are stored or corrosive gases exist;

### Product Packaging List

| Name                     | Image Quantity |
|--------------------------|----------------|
| Main Unit                | 1              |
| Filter                   | 2              |
| Dermabrasion Handle      | 1 box          |
| Transparent Suction Tube | 1              |
| Power Cord               | 1              |

### Product Specifications

- Input Voltage: AC 220V-240V
- Product Frequency: 50/60HZ
- Product Power: 20W
- Packaging Dimensions: 22x22x20CM
- Gross Weight: 2.4KG
- Net Weight: 2.15KG



1. Connect the device to a 220V power supply and press the power button. The device will emit a "beep" sound, indicating that it is powered on and enters standby mode;
2. Insert the suction tube into the filter nozzle, connect the suction pen to the suction tube, and insert the selected dermabrasion head. Press  the small hole on the handle to use it normally;
3. Adjust the suction power from low to high using the  knob.

14. This instrument shall not be placed or used near radios, televisions, copiers, fax machines, etc.;
15. Do not use it together with other instruments (such as microwaves) to avoid mutual interference between instruments and cause misoperation;
16. If an abnormal situation occurs to the instrument or the user during operation, the operator shall immediately move the operation handle away from the user, press the stop switch, and stop using the instrument;
17. The instrument shall be correctly used by skilled operators. Personnel other than operators are not allowed to use the instrument;
18. After the operation, place the control handle back on the control handle holder;
19. After using the operation handle, clean the probe thoroughly for next use.

### Precautions for Troubleshooting, Storage, and Inspection

1. Do not use the instrument if it malfunctions. In case of smoke, burning smell, or inability to operate, continued use may cause fire or electric shock. Immediately press the stop button, turn off the power, unplug the power plug, and contact the company or distributor immediately;
2. Personnel other than professional technicians are strictly prohibited from repairing or modifying the instrument. The interior of the instrument must never be opened;
3. The instrument and its components must be inspected regularly;
4. Do not use cleaning tools that may damage the instrument's protective coating, such as brushes, brooms, etc.;
5. If the instrument is not used for more than 2 weeks, it must be carefully inspected before reuse.

## Contraindicated Groups

### (Do Not Use If You Belong to the Following Groups)

- Individuals with implanted cardiac pacemakers, defibrillators, or metal implants (such as steel plates, bone nails) in the body;
- Pregnant women, lactating women, and women during menstruation;
- Patients with epilepsy, hypertension, severe heart disease, diabetes, tumors or cancer, or thyroid diseases;
- Individuals with open wounds, skin damage, severe inflammation, eczema, dermatitis on the skin, or those with facial telangiectasia (spider veins) or rosacea;
- Those who have received botulinum toxin injections, filler injections, laser treatments, chemical peels, or plastic surgery recently (usually within 6 months);

### **Individuals with dry or sensitive skin, or those with damaged, inflamed (such as severe acne, eczema), or sunburned skin.**

#### Working Principle and Operation of the Diamond Dermabrasion Machine

##### Working Principle

The instrument uses a probe inlaid with fine diamond particles to perform physical exfoliation on the skin, gently removing aging cuticles. Combined with vacuum negative pressure technology, it suctions away the exfoliated cuticles and dirt in the pores.

##### Operation Process

1. Thoroughly cleanse the face;
2. Apply a warm towel to the face for 1-2 minutes to help soften cuticles and open pores;
3. Use the supporting spray or gel as a medium to reduce friction between the instrument and the skin.

#### Zonal Operation

- T-Zone (Forehead, Nose, Chin): Usually has more oil and blackheads, suitable for using suction mode or low-grade dermabrasion. On blackhead-prone areas such as the nasal alae, move slowly, but do not repeat suction on the same area more than 3 times;
- U-Zone (Cheeks): The skin here is thinner, suitable for lifting and gentle dermabrasion. It is recommended to slide from bottom to top, against the direction of wrinkles.

#### Post-Use Repair

After use, immediately soothe and repair the skin. You can apply a repair mask containing ingredients such as Vitamin B5 and Ceramide. Afterwards, enhance moisturization and strictly protect against the sun, as the newly formed cuticle is more sensitive to ultraviolet rays.

#### Usage Frequency

For home use, the diamond dermabrasion machine should not be used more than 2-3 times a week. Allow the skin sufficient time for self-repair to avoid damage to the skin barrier caused by over-cleansing.

#### Important Safety and Maintenance Instructions

- Force and Duration: During operation, be sure to be gentle and slow. Rely on the instrument's own functions; do not press hard or stay in the same position for a long time. The recommended single full-face use time is controlled within 10-15 minutes;
- Instrument Maintenance: After use, clean the probe and filter in a timely manner as required by the manual. The diamond probe is for personal use only; it is recommended to replace it regularly to avoid bacterial growth. Store the instrument in a dry and cool place.