your wellness company


EN

Instructions Manual

Beauty instrument

F807A.N

VACCUM & SPRAY

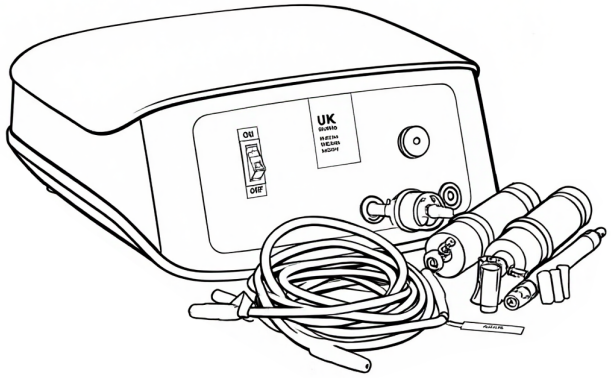

Imported by
WEELKO BARCELONA SL

Pol. Ind. La Valldan
Camí de Garreta, 92
08600 Berga
Barcelona - Spain
T. +34 93 824 90 20
E. info@weelko.com

VERSION 1.0

weelko your wellness company

Thank you sincerely for choosing our product.. This manual is compiled to help
users operate this instrument ef_fectively. We hope our company's products and
services will bring you a distinguished experience.

This user manual contains all general information about the instrument. It serves
as the starting point for new users to begin using the instrument, as it provides a
general introduction to the system, explains how it works, and details the instru-
ment's installation and various operation methods.

Usage Precautions

1. Before turning on the machine, be sure to read the user manual carefully. Operators must be familiar with the user manual and operation methods, and the instrument must be managed by designated person;
2. When replacing the fuse, always turn off the power switch to prevent f_ire and electric shock;
3. Use a power supply of 220V±22V, 50HZ±1HZ and an ef_fectively grounded socket;
4. Using a damaged or modif_ied power cord may cause f_ire and electric shock;
5. If the instrument is loose, falling, or damaged, please contact the company in time;
6. Do not use the instrument in unstable environments such as tilting, vibrating, or impacting conditions;
7. No water or debris is allowed to enter the interior of the instrument. In case of foreign objects entering, please contact the company promptly;
8. When moving the instrument, please unplug the power cord first;
9. When moving the instrument, please unplug the power cord first;
10. The instrument should be placed on a stable surface, with at least 5cm of space left around it to ensure air circulation and facilitate heat dissipation of the instrument;
11. Do not cover the instrument with anything to prevent blockage of the air vents;
12. Avoid using the instrument in the following environments: direct sunlight; high humidity environments; near water; vibrating or tilting places; strong magnetic f_ield environments; places where electromagnetic waves or impulse voltages occur; places where chemicals are stored or corrosive gases exist;

Product Packaging List

NameImage	Quantity
Main Unit1	
Filter1	
Suction Pen (with 3 suction heads)	1 box
Spray Bottle2	
Power Cord1	

Product Specif_ications

Input Voltage: AC 220V-240V
Product Frequency: 50/60HZ
Product Power: 60W
Packaging Dimensions: 34x34x22CM
Gross Weight: 2.75KG
Net Weight: 2.25KG

14. This instrument shall not be placed or used near radios, televisions, copiers,
15. fax machines, etc.;
    Do not use it together with other instruments (such as microwaves) to
16. avoid
    mutual interference between instruments and cause misoperation;
    If an abnormal situation occurs to the instrument or the user during
17. operation,
    the operator shall immediately move the operation handle away from the
18. user, press the stop switch, and stop using the instrument;
    The instrument shall be correctly used by skilled operators. Personnel other
19. than operators are not allowed to use the instrument;
    After the operation, place the control handle back on the control handle holder;
    After using the operation handle, clean the probe thoroughly for next use.

Precautions for Troubleshooting, Storage, and Inspection

1. Do not use the instrument if it malfunctions. In case of smoke, burning smell,
   or inability to operate, continued use may cause f_ire or electric shock.
   Immediately press the stop button, turn of_f the power, unplug the power plug, and contact the professional technician immediately;
   Please do not disassemble the instrument by yourself. It is strictly prohibited
2. from
   repairing or modifying the instrument. The interior of the instrument must The instrument and its components must be inspected
   Do not use sharp tools that may damage the instrument's regularly;
3. protective
   coating, such as brushes, brooms, etc.;
   If the instrument is not used for more than 2 weeks, it must be
4. carefully
   inspected before reuse.

Contraindicated Groups (Do Not Use If You Belong to the Following Groups)

Individuals with implanted cardiac pacemakers, def_ibrillators, or metal implants (such as steel plates, bone nails) in the body;

Pregnant women, lactating women, and women during menstruation;

Patients with epilepsy, hypertension, severe heart disease, diabetes, tumors or cancer, or with infectious diseases;

Individuals with open wounds, skin damage, severe inf_lammation, eczema, dermatitis on the skin, or those with facial telangiectasia (spider veins) or rosacea;


Those who have received botulinum toxin injections, f_iller injections, laser treatments, chemical peels, or plastic surgery recently (usually within 6 months); or

Individuals with skin breakdown, inf_lammation, severe acne, eczema, or rosacea.


Please pay attention to the following:

1. Frequency Control: Do not use it daily. It is recommended to use the deep cleansing mode 1-2 times a week; the specif_ic frequency depends on individu-
   al skin type.
2. Intensity and Duration: Prioritize lower intensity and shorter duration. Do not
   use high intensity for prolonged periods.
3. Post-Treatment Recovery Period: After facial plastic surgery, laser treatment,
   chemical peeling, and other projects, wait for the skin to fully recover (usually
   at least 2 months is recommended) before using the instrument;
4. Special Reminder for Pregnant Women: Due to unstable hormone levels, preg-
   nant or lactating women are not recommended to use the instrument.

Instrument Introduction

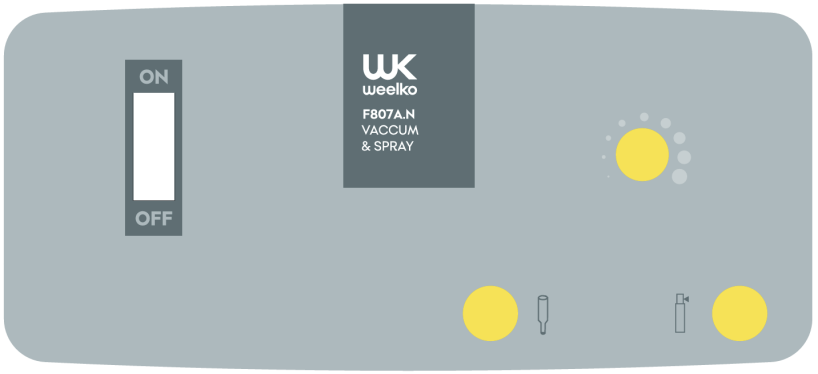

Instrument Operation Steps

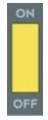

1. Connect the device to a 220V power supply and press the power button .
   The device will emit a "beep" sound, indicating that it is powered on and enters standby mode;
2. For Suction Mode: Insert the suction tube into the f_ilter nozzle, connect the suction pen to the suction tube, and insert the selected suction head.
3. For Spray Mode: Insert the suction tube into the spray bottle connector and press the small hole on the handle to use it normally;
   Press the small hole of the air vent above the spray bottle to atomize the product;
4. Suction Adjustment: Adjust the suction power from low to high using the knob.