your wellness company


EN

Instructions Manual

Beauty instrument

F806.N

HIGH FREQUENCY

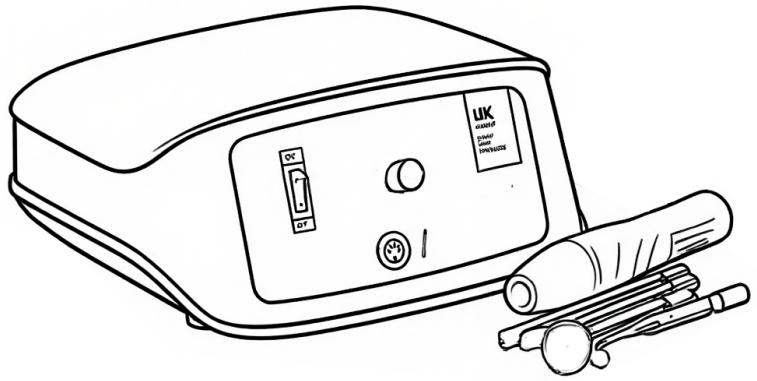

Imported by
WEELKO BARCELONA SL

Pol. Ind. La Valldan
Camí de Garreta, 92
08600 Berga
Barcelona - Spain
T. +34 93 824 90 20
E. info@weelko.com

VERSION 1.0

weelko your wellness company

Thank you sincerely for choosing our product.. This manual is compiled to help
users operate this instrument ef_fectively. We hope our company's products and
services will bring you a distinguished experience.

This user manual contains all general information about the instrument. It serves
as the starting point for new users to begin using the instrument, as it provides a
general introduction to the system, explains how it works, and details the instru-
ments installation and various operation methods.

Usage Precautions

1. Before turning on the machine, be sure to read the user manual carefully. Operators must be familiar with the user manual and operation methods, and the instrument must be managed by designated person;
2. Do not use a power supply that does not meet the specified voltage;
3. When replacing the fuse, always turn off the power switch to prevent f_ire and electric shock;
4. Use a power supply of 220V±22V, 50HZ±1HZ and an ef_fectively grounded socket;
5. Using a damaged or modif_ied power cord may cause f_ire and electric shock hazards;
6. If the instrument is loose, falling, or damaged, please contact the company in time;
7. Do not use the instrument in unstable environments such as tilting, vibrating, or impacting conditions;
8. No water or debris is allowed to enter the interior of the instrument. In case of foreign objects entering, please contact the company promptly;
9. When moving the instrument, please power off and unplug the power cord f_irst;
10. When moving the instrument, try to avoid vibration;
11. The instrument should be placed on a stable surface, with at least 5cm of space left around it to ensure air circulation and facilitate heat dissipation of the instrument;
12. Do not cover the instrument with anything to prevent blockage of the air vents;
13. Avoid using the instrument in the following environments: direct sunlight; high humidity environments; near water; vibrating or tilting places; strong magnetic f_ield environments; places where electromagnetic waves or impulse voltages occur; places where chemicals are stored or corrosive gases exist;

Product Packaging List

Name	Image	Quantity
main machine		1
High frequency handle		1
Electrotherapy tube (tongue shaped)		1
Electrotherapy tube (straight shaped)		1
Electrotherapy tube (curved shaped)		1
Electrotherapy tube (mushroom shaped)		1
power cord		1

Product Specif_ications

Input Voltage: AC 220V-240V
Product Frequency: 50/60HZ
Product Power: 60W
Packaging Dimensions: 34x34x22CM
Gross Weight: 2.75KG
Net Weight: 2.25KG

14. This instrument shall not be placed or used near radios, televisions, copiers,
15. fax machines, etc.;
    Do not use it together with other instruments (such as microwaves) to
16. avoid
    mutual interference between instruments and cause misoperation;
    If an abnormal situation occurs to the instrument or the user during
17. operation,
    the operator shall immediately move the operation handle away from the
18. user, press the stop switch, and stop using the instrument;
    The instrument shall be correctly used by skilled operators. Personnel other
19. than operators are not allowed to use the instrument;
    After the operation, place the control handle back on the control handle holder;
    After using the operation handle, clean the probe thoroughly for next use.

Precautions for Troubleshooting, Storage, and Inspection

1. Do not use the instrument if it malfunctions. In case of smoke, burning smell,
   or inability to operate, continued use may cause f_ire or electric shock.
   Immediately press the stop button, turn of_f the power, unplug the power plug, and contact the professional technician immediately;
2. Please do not disassemble, repair or modify the instrument by yourself. Do not use prohibited from
   repairing or modifying the instrument. The interior of the instrument must The instrument and its components must be inspected
   Do not use sharp tools that may damage the instrument's regularly;
3. protective
   coating, such as brushes, brooms, etc.;
4. If the instrument is not used for more than 2 weeks, it must be
5. carefully
   inspected before reuse.

Contraindicated Groups (Do Not Use If You Belong to the Following Groups)

Individuals with implanted cardiac pacemakers, def_ibrillators, or metal implants (such as steel plates, bone nails) in the body;

Pregnant women during menstruation;

Patients with epilepsy, hypertension, severe heart disease, diabetes, tumors or cancer, or with infectious diseases;

Individuals with open wounds, skin damage, severe inf_lammation, eczema, dermatitis on the skin, or those with facial telangiectasia (spider veins) or rosacea;

Those who have received botulinum toxin injections, f_iller injections, laser treatments, chemical peels, or plastic surgery recently (usually within 6 months).

Operation Process

1. Pre-Use Preparation

Inspect the Instrument: Ensure the instrument is dry and intact, and the glass electrode tube is free of cracks.

Clean the Skin: Thoroughly cleanse the treatment area to remove makeup and impurities, ensure the treatment area is dry and free of oil.

Connect the Electrodes: First select and insert the appropriate glass electrode, then apply a sufficient amount of special conductive gel or essence on the skin. Dry operation is strictly prohibited to provide conductivity and avoid burns.

2. Operation Process

Power On and Adjustment: Turn on the power and be sure to start from the lowest intensity gear. Usually, rotating clockwise increases the intensity, and rotating counterclockwise decreases it;


Control Duration: Do not stay in the same area for more than 10 seconds to avoid damage;

Avoid Dangerous Areas: Strictly prohibit use on eyeballs, throat, open wounds, severe acne or infected areas.

3. Post-Use Operation

Shutdown Sequence: First adjust the intensity knob back to the lowest gear or position, then turn off the power, and unplugging the power, carefully remove the glass electrode tube from the instrument handle, and clean the glass electrode tube with a soft damp cloth, and store them in a cool and dry place after drying.

Instrument Introduction

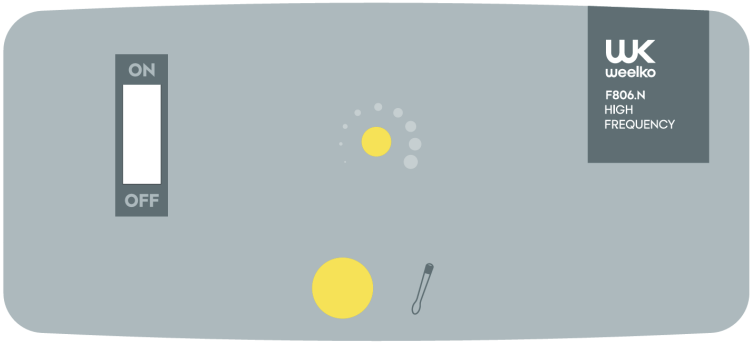

Instrument Operation Steps

1. Connect the device to a 220V power supply and press the 

power button.
The device will emit a "beep" sound, indicating that it is powered on and enters standby mode;
2.  Insert the high-frequency handle according to the accessory mark, and install the selected electrode tube into the handle;
3. Adjust the intensity adjustment knob  to the high or low level as needed;
4. After setting, start the instrument to begin treatment. During use, follow the operation techniques and duration requirements in the "Operation Process" section;
5. After the treatment, follow the "Post-Use Operation" steps to shut down, disassemble, clean and store.