EN

SILVERFOX

2296 Series Facial Chair

Welcome to use our product.
The manufacturer reserves the right to improve design and modify the product without prior notice.
Version A0

1. About This Manual

This manual is for reference only when operating this specific model. If the manual is used for other purposes, we will not undertake any resulting consequences or responsibilities.
All information of this manual shall be protected by copyright and shall not be copied and translated in any part of this specification without the written consent of our company SILVERFOX.

2. Warning

Electrical Hazards

- The plug must be plugged into an outlet that is properly installed and with protective earth, to prevent the risk of fire and electric shock.
- The ground wire in the power cord serves only as a functional ground.
- Do not place the device in a location where it is difficult to disconnect the power supply.
- The power supply must be cut off when replacing electrical components or cleaning the equipment.
- Do not place the device in humid places such as bathrooms to avoid electric shock or other malfunctions; keep the hand remote controller and control box away from water.
- If the plug is damaged, please contact the dealer, service center or relevant personnel with relevant skills and experience for repair to avoid danger.

Mechanical Risk

- Do not stand on this equipment.
- Do not allow children to operate or play with and place with the equipment or place it within their reach.
- The maximum load capacity of this chair is 200 kg. Exceeding the load capacity may cause safety risks.

Environmental Controls

- Please do not place and operate the equipment in places where is painting/coating or supplying oxygen.
- Electric shock protection: Class I
- For indoor use only

Operational Protocol

- After usage, please stay the upholstery flat to 180°.
- This equipment should be operated by trained personnel who are familiar with its structure, performance, and operating instructions, and who should perform regular maintenance.
- This equipment shall not be used by persons with reduced physical, sensory, or mental capabilities without adequate experience or knowledge, unless they are under continuous supervision or have received instructions regarding safe operation by a qualified person responsible for their safety."
 WARNING: During equipment operation, verify that the range of motion is clear of obstacles to prevent crushing or entanglement hazards.

Disposal & Modification

- Any modifications to this device are not permitted.
- The equipment and accessories must be disposed of according to local regulations after their service life ends.

3. Precautions


Maintenance Protocol

- Do not repair or maintain this device during treatment operations.
- If the user and/or patient becomes aware of any serious incident associated with this equipment, it should be reported to the regulatory authorities in the user/patient's country and to the manufacturer.
- Regularly tighten bolts and nuts to prevent malfunctions caused by loose components, slippage, vibration, and falling.
- During maintenance, operators must not disassemble any components of the chair.
- Do not operate this product if it malfunctions. If any abnormal phenomenon or sound is detected, immediately stop operation and inspect it. Operation can only be resumed after the malfunction is corrected.
- Maintenance should be performed by trained professionals.

- After each treatment, use a neutral solution to remove dust and stains from the leather.
- Replacements must be performed by trained professionals. Our company can provide circuit diagrams, troubleshooting guides, or other helpful maintenance materials for professional reference.


Operational Protocol

- Please read this manual carefully before installation, operation, and adjustment, and keep it in a safe place.
- Observe safety operating procedures, do not overload the table, ensure that all safety devices are complete and reliable, and promptly eliminate any unsafe conditions.
- Operate the product in strict accordance with the warning labels attached to the product. To control the entire chair's lift, the button must be pressed continuously until the chair is in the desired position.
- Avoid single-point load operations. (distribute weight evenly)
- Motor duty cycle: ≤2minutes ON / ≥18minutes OFF. (the motor is designed to operate for a maximum of 2 minutes, followed by a minimum of 18 minutes of rest)
- When adjusting each component, keep your hands, feet, and head away from the arrow points that shown in below picture.

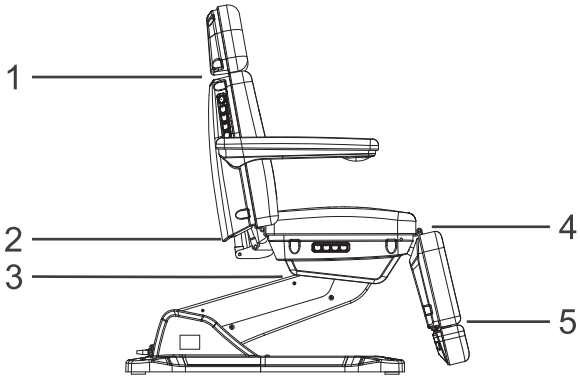4. Installation

- Unpack the box and inspect the product to ensure it is in good condition and that all accessories and spare parts are present. If you have any questions, please contact your dealer.
- Place the product in a dry, clean, treatment-friendly area of the room, leaving ample space around the device for operation. Ensure the floor is level. If the floor is uneven, rotate the rubber feet to compensate for any unevenness to ensure stability.
- Note: No additional screws are required to secure the table to the floor.
- Connect the product plug to the power source properly (as shown in the image).

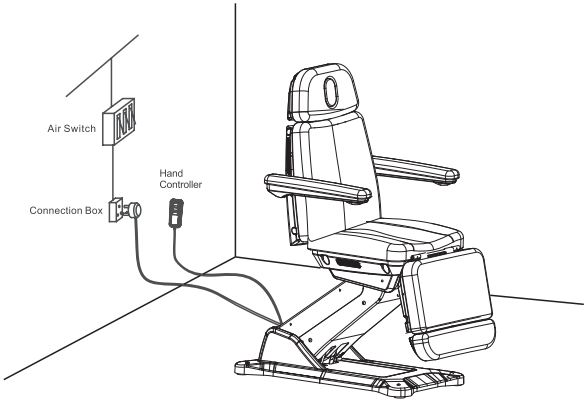5. Structure of products


This device consists of a control box, motor, hardware, plastic components, and leather.

Function description


Electrically operated beauty bed, powered by an internal motorized system for precise and adjustable positioning. Designed to provide stable support for clients during various skincare and beauty treatments.

Circuit System Diagram

2296

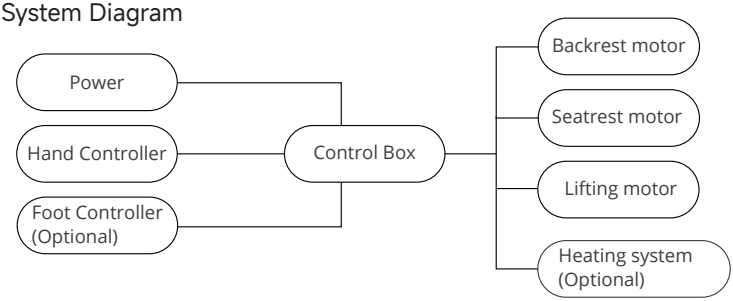
2296B
2296B PRO

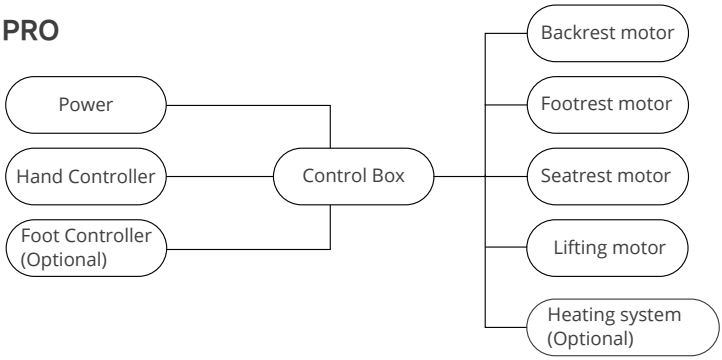The control box


Depending on the model, this device includes backrest angle adjustment, seat cushion angle adjustment, and lift adjustment. It is an overall integral structure. (As shown in picture)

2296

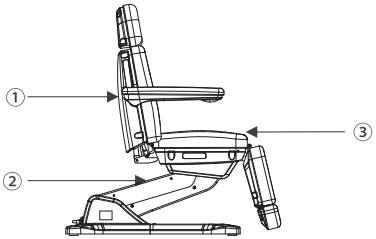① angle adjustment of the backrest
② lifting adjustment
③ angle adjustment of the seat
④ angle adjustment of the footrest (Manual adjustment)

2296B / 2296B PRO

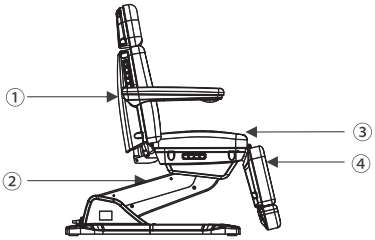
① angle adjustment of the backrest
② lifting adjustment
③ angle adjustment of the seat
④ angle adjustment of the footrest


6. Operation Instruction

2296


1.Standard Configuration
(Hand Controller)

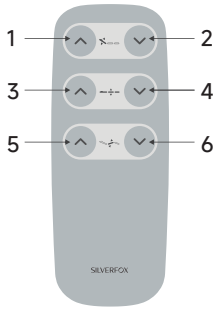2.Optional Configuration
(Hand Controller)

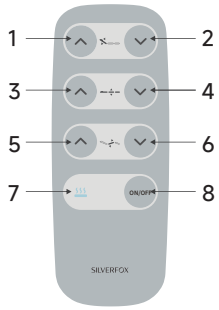
3. Optional Configuration
(Foot Controller)

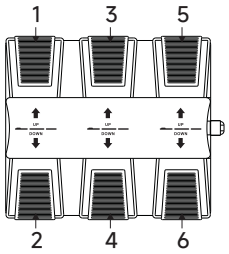
No.	Description	Operation
1/2	Backrest UP or DOWN	Press and hold the appropriate button to move the chair Movement stops when… Button is released Maximum / minimum position is reached.
3/4	Chair Lift UP or DOWN	
5/6	Tilt UP or DOWN	
7	Indicator light of heating function	Pressing the button, the light is on, the heating function works; pressing the button again, the light is off, the heating function stops working.
8	Heating button(Optional)	

2296B


1.Standard Configuration
(Hand Controller)

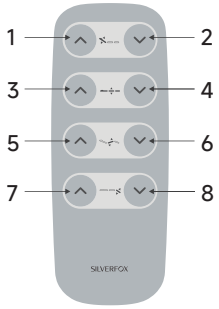2.Optional Configuration
(Hand Controller)

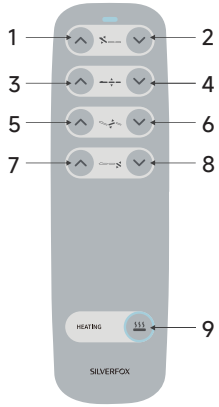
3. Optional Configuration
(Foot Controller)

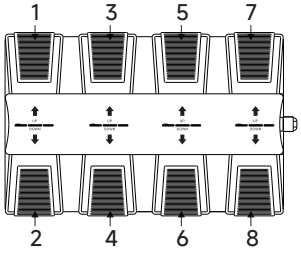
No.	Description	Operation
1/2	Backrest UP or DOWN	Press and hold the appropriate button to move the chair Movement stops when… Button is released Maximum / minimum position is reached.
3/4	Chair Lift UP or DOWN	
5/6	Tilt UP or DOWN	
7/8	Foot UP or DOWN	
9	Heating button (Optional)	Pressing the button, the light is on, the heating function works; pressing the button again, the light is off, the heating function stops working.

2296B PRO


1.Standard Configuration
(Hand Controller)


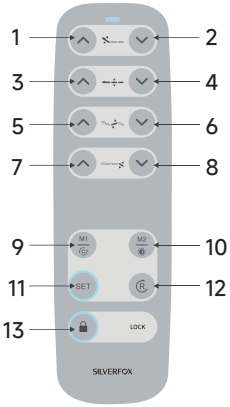
2.Optional Configuration
(Hand Controller)

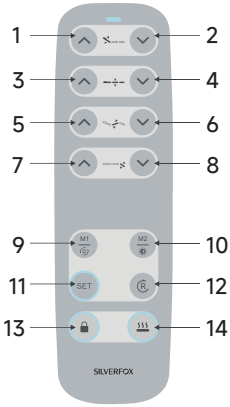3. Optional Configuration
(Foot Controller)

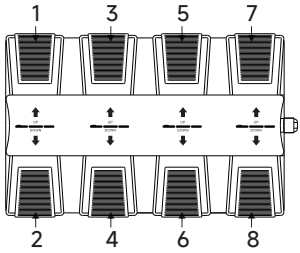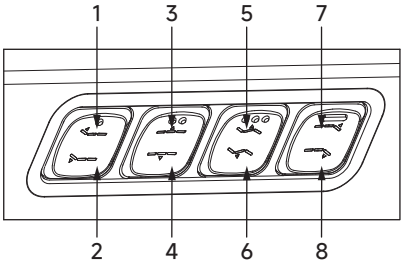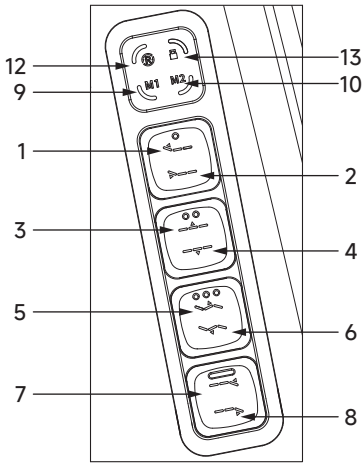No.	Description	Operation
1/2	Backrest UP or DOWN	Press and hold the appropriate button to move the chair Movement stops when… Button is released Maximum / minimum position is reached.
3/4	Chair Lift UP or DOWN	
5/6	Tilt UP or DOWN	
7/8	Foot UP or DOWN	
9/10 Upper Section	 Memory Function Programming the chair positions	The Programming feature allows you to easily recall two frequently used chair positions with M1, M2 buttons. -Adjust the chair to the desired position -Press and hold the "M1" or "M2" button (on the backrest button or the hand remote control) for 3 seconds. You will hear a beep and the position will be set to "M1" or "M2". To recall a programmed position… -Press the desire position button “M1” or “M2” Movement stops when… -Programmed position is reached -Any button is pressed (on backrest, seatrest, foot pedal or hand remote controller)
9/10 Lower Section 11	 Ambient Light	There are 14 lighting modes and 12 brightness levels. Briefly press the “SET” button. The “SET” indicator will flash, indicating entry into the light strip control mode. Briefly press the “SET” button again, or leave the light strip buttons inactive for more than 10 seconds, to exit the mode. When the light strip is OFF, briefly press “M1” to turn it ON. When the light strip is ON: -Long press “M1” to turn the light strip OFF. -Briefly press “M1” to switch between lighting modes: 1. Solid Color Modes: White, Red, Green, Blue, Red + Green, Red + Blue 2. Single Color Fade Modes: White Fade, Red Fade, Green Fade, Blue Fade 3. Single Color Flash Modes: White Flash, Red Flash, Green Flash, Blue Flash 4. Lights OFF “M2” is used for brightness adjustment of the light strip, with a total of 12 levels. Brightness can be adjusted from dim to bright or from bright to dim. When the minimum or maximum level is reached, the light strip will flash twice as an indication, then reverse the adjustment direction. Both short press and long press operations are supported.

12	Reset	At any position, press “Reset” button (on the backrest, seatrest, or hand remote controller) … -The chair automatically moves to the entrance position so that your patient may get on or off the chair. Backrest moves to the upright position Chair tilt moves to the flat position Chair height moves to the minimum foot moves down to the maximum Movement stops when… -Entrance position is reached -Any button is pressed (on backrest, seatrest, foot pedal or hand remote controller) *Note: To avoid damaging the chair, the reset function will be disabled when the footrest is extended. Before using the reset function, ensure that the foot extension is retracted and firmly pressed against the footrest. At this point, the safety switch of the extension device is activated, and the reset function will take effect.
13	Lock	Press and hold for 3 seconds to activate the lock function. When the lock mode is activated, all other buttons are disabled and the motors will not operate. Press and hold for 3 seconds again to deactivate the lock function.
14	Heating button (Optional)	Pressing the button, the light is on, the heating function works; pressing the button again, the light is off, the heating function stops working.


Notes:
1.The heating function is only available for the models with heating function.
2.Different hand controllers should be chosen according to different configuration.


7. Technical Parameters


Power	100-120V～ or 220-240V～ ,50/60Hz
Maximum Power	300W
Environment Temperature	0°C ~ + 40°C
Relative Humidity	30% ~ 90%
Lifting Rated Load	200kg
Package Size	2296 :1500 × 620 × 615 mm 2296B / 2296B PRO:1370 × 620 × 615 mm

2296 / 2296B / 2296B PRO

Unit: mm

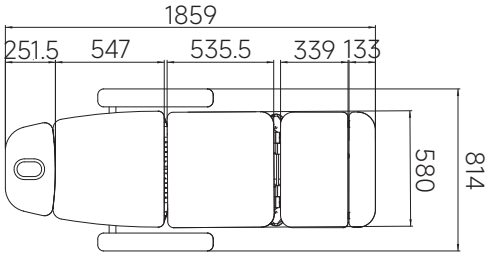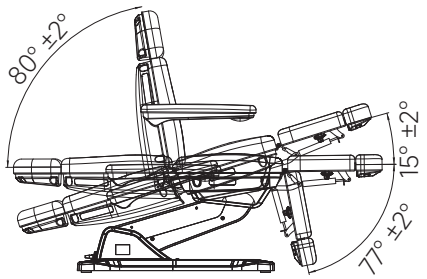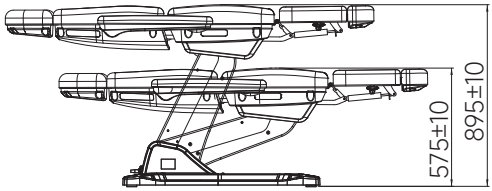8. Maintenance

1. After every treatment, please only use neutral solution to clear the dust, stains on leather, metal part or other parts. Disinfecting solution like alcohol, acid or alkaline solution are forbidden.This is to prevent the lifetime of the components from being shortened.
2. Please lubricate the moving parts regularly (once a year).
3. When the headrest is stuck, it can be improved by applying lubricating oil.

9. Troubleshooting

The Fault Phenomenon	Reason	Solution
motor can't work normally	1. The power cord is unconnected or not plugged well.	1. Connect the power cord or plug it well.
	2. Poor quality of the hand controller or foot controller.	2. Replace the hand controller or foot controller.
	3. Poor quality of the control box.	3. Please contact maintenance center or supplier for troubleshooting.
	4. Repture of motor tail.	
	5. Motor was damaged.	
Abnormal Sound	1. The rotating parts are short of lubricating oil.	1. Add lubricating oil.
	2. Rotating positions (such as hinge, POM sets, gas lift) are short of lubricating oil.	2. Add lubricating oil.
	3. Motor quality problems.	3. Replace the motor.

If the problem persists after following the above steps
please contact maintenance center or supplier to solve it.

10. Disposal

Contact its local authorities to determine the proper method of disposal.

11. Symbols explanation

	Refer to instruction manual		Manufacturer
	WEEE		CE
	Date of manufacturer		Authorized representative in the European Community
	Warning		Series NO.
	Keep dry from rain rom rain		Fragile, handle with care
	Keep away from sunlight		Stacking limit by mass of 250kg
	Upward		No tumbling
	Be careful		No clasping

CE

The machinery fulfils the following relevant provisions:

Directive 2006/42/EC

Directive 2014/30/EU

Directive 2014/35/EU

Directive2011/65/EU

SILVERFOX CORPORATION LIMITED
Adresse: NO.18, 1st TongLe Road, TangXia Town, PengJiang District, JiangMen City, GuangDong Province, P.R. China.

According to Waste of Electrical and Electronic Equipment (WEEE) directive, WEEE should be separately collected and treated. If at any time in future you need to dispose of this product please do NOT dispose of this product with household waste. Please send this product to WEEE collecting points where available.

▲ Product components may be subject to changes without prior notice. Actual appearance may vary from images.

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