F-366

Beauty Instrument

User Manual

Welcome to use our products. A0
Please read the user manual carefully before using the products.
Our company reserves the right to improve design and modify
the products without prior notice.

Catalog

About This Manual ---- 1

Safety Guidance ---- 1

Components ---- 2

Product Parameters ---- 3

Panel and Function Introduction ---- 4

Function Introduction ---- 5/8

Troubleshooting ---- 9

About This Manual

This manual shall only be used for the reference of operation of this particular model. If the manual is used for other purposes, SILVERFOX will not undertake any according consequences and responsibility. All information of this manual shall be protected by copyright and shall not be copied and translated in any part of this specification without the written consent of our company SILVERFOX.

Safety Guidance

To ensure your correct and safe use of this product, please read all the instructions before use! Please follow all safety precautions!


	**Warning!**	**Indicates a potential dangerous situation that, if not avoided, could result in serious loss of life.** • This product is only for beauty practitioners with professional qualifications. • Do not use or store the product in a damp place such as bathroom to avoid malfunction and danger. • Do not use this product near flammable items (e.g., alcohol, thinner, gasoline, nail polish remover, sprays, etc.). • Do not use this product under the following conditions: when the power cord or plug is damaged or gets hot; when the plug in the power socket comes loose. • Do not use wet hands to plug and unplug power cord connectors and use this product. • Please make sure that the working voltage of this product conforms to its rated voltage. • Before use, make sure that the power plug and the power cord connector are fully plugged in, so as not to cause danger. • When not using this product, please turn off the power and pull out the plug. • Please regularly clean the dust accumulated on the power plug and power cord connector to avoid failure of insulation caused by moisture, dust and other factors. • Be sure to use the special power cord assembly for this product. Do not use the power cord assembly on any other appliance. • Do not wrap the power cord around the product, otherwise it may distort and malfunction the power cord. • Keep out of reach of children. Do not allow children to use this product. • Please do not modify, disassemble or repair this product by yourself. If you need to repair it, please contact the manufacturer, agent or professional maintenance personnel to check and repair it, so as not to cause danger. • Please contact distributor, service center or related personal with appropriate skill and experience for repairing if the plug is damaged, in case of any danger. • Please be sure to use this product according to the purpose specified in this manual. Do not use accessories not recommended by the manufacturer.
	**Attention!**	**Represents a potentially dangerous situation which, if not avoided, can result in general injury to person.** • Operators should pay special attention to the safety of personnel and equipment when using. • If the product is abnormal, please stop using it immediately and contact the manufacturer or dealer for relevant treatment. • This product and related accessories should be handled lightly. • Please properly handle the packaging bag of this product to avoid accidents such as asphyxiation of infants, and keep it out of reach of children.

◆ F-366




1. Probe for face
2. Probe for eye
3. High-Frequency

Rear Panel

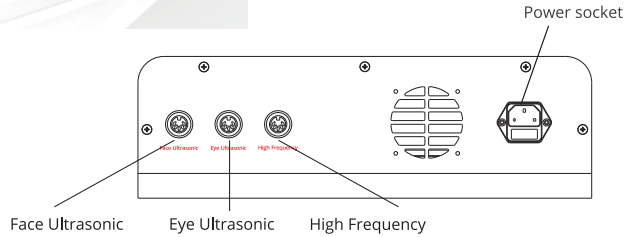

◆ Front Panel

1. Face Ultrasonic
2.Eye Ultrasonic
3.High Frequency

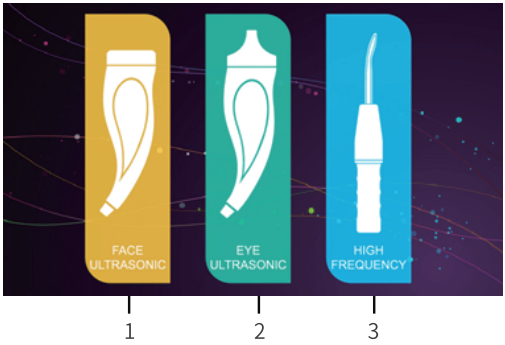

◆ Operation Method

Operation procedure

① Click to choose the required function and it will enter the interface of the corresponding function.


② All working interfaces are with the button to return to the main interface.

Product Parameters

Product model	F-366
Rated voltage / Frequency	AC 100-120V/220-240V,50/60Hz
Power	22W
Net weight	6.1kg
Product size	37×31.5×16.5cm
Package size	56.8×52.5×23.5cm
Package contents	Face Piece for Ultrasonic× 1
	Eye Piece for Ultrasonic× 1
	High-Frequency× 1
	Main Unit× 1
	Accessories Shelf× 1
	Power Cord× 1

Introduction

1. Face Ultrasonic
2.Eye Ultrasonic
3.High Frequency

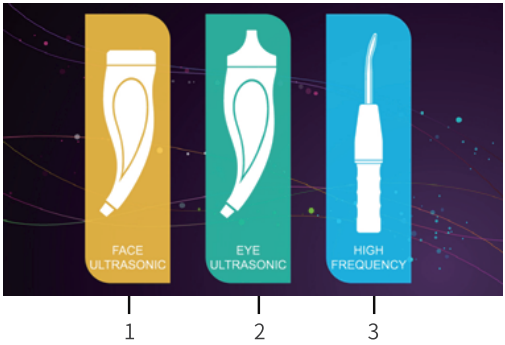

Operation Method

1.Click to choose the required function and it will enter the interface of the corresponding function.
2.All working interfaces are with the button to return to the main interface.

◆ Operation Cautions

⚠ Warning

- Do not touch the surface of the accessory probe and the output interface when the machine is powered on.
- Do not point the probe directly at the eyeball or touch the eye.
- It is forbidden for the following groups:
  1.Allergic sufferer
  2.People with skin burns, scalds, sunburns, inflammation or other wounds
- It is forbidden to start to use the probe without cosmetic operation to prevent the risk of high temperature scald and prevent to reduce the service life of the probe due to overheating.
- Do not remove the probe and handle by yourself. If turn the probe and handle, the probe may be damaged.

⚠ Attention

- The probe should be kept clean at all times and disinfected with 75% alcohol after each use.
- Avoid prolonged contact of the probe with the skin in the same position. If the operation time is too long, it will not only fatigue the skin and affect the absorption, but also may cause danger due to the mechanical vibration generated by the output to increase the temperature of the probe.


Recommended using time

Adjust the time according to the specific situation of the nursing site and skin, generally the whole face operation time is 10-15 minutes.

Function Introduction



Face Ultrasonic probe (Cosmetic use only)

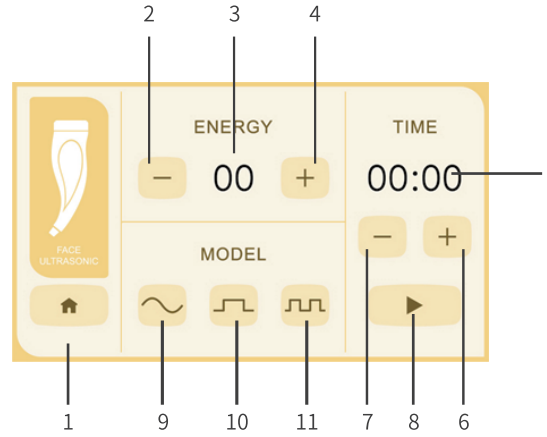

1. Return to main interface
2. Reducing Energy Level
   Energy level range:1~10
3. Energy Level Display
4. Increasing Energy Level
   Energy level range:1~10
5. Working time display
   The working time enters the countdown when operating. The machine will stop automatically when the time shows zero.
6. Increasing working time
   Working time range: 1~30 minutes.
7. Reducing working time
   Working time range: 1~30 minutes.
8. Start/Pause
9. Continuous Wave Mode
   Output energy is continuous.
10. Intermittent wave mode
    Output energy is intermittent.
11. Mixed Wave Mode
    Switching of continuous wave mode and mixed wave mode.

Operation Method


1) Setting energy output level and working time. The maximum time can be set to 30min.
2) Select energy output mode
3) Take out the ultrasonic handle and turn on the “start” button
4) The ultrasonic stops automatically when scheduled time ends. Clean the metal probe and put it back to the accessories shelf.

- After the nursing operation, the probe shall be disinfected with 75% alcohol and properly dry stored.

Note: Do not wet the metal part of the probe.

Eye Ultrasonic probe (Cosmetic use only)



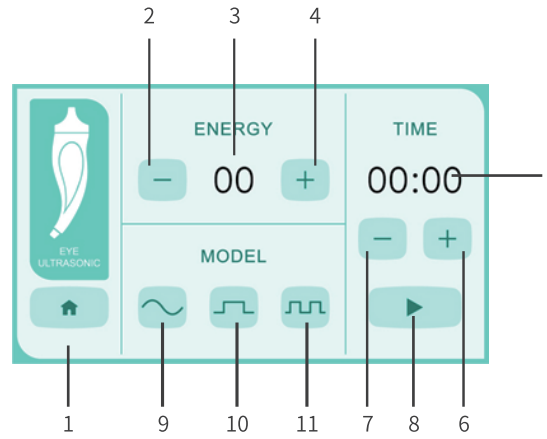

1. Return to main interface
2. Reducing Energy Level
   Energy level range:1~10
3. Energy Level Display
4. Increasing Energy Level
   Energy level range:1~10
5. Working time display
   The working time enters the countdown when operating. The machine will stop automatically when the time shows zero.
6. Increasing working time
   Working time range: 1~30 minutes.
7. Reducing working time
   Working time range: 1~30 minutes.
8. Start/Pause
9. Continuous Wave Mode
   Output energy is continuous.
10. Intermittent wave mode
    Output energy is intermittent.
11. Mixed Wave Mode
    Switching of continuous wave mode and mixed wave mode.

Operation Method


1) Setting energy output level and working time. The maximum time can be set to 30min.
2) Select energy output mode
3) Take out the ultrasonic handle and turn on the “start” button
4) The ultrasonic stops automatically when scheduled time ends. Clean the metal probe and put it back to the accessories shelf.

- After the nursing operation, the probe shall be disinfected with 75% alcohol and properly dry stored.

Note: Do not wet the metal part of the probe.



High frequency electrotherapy

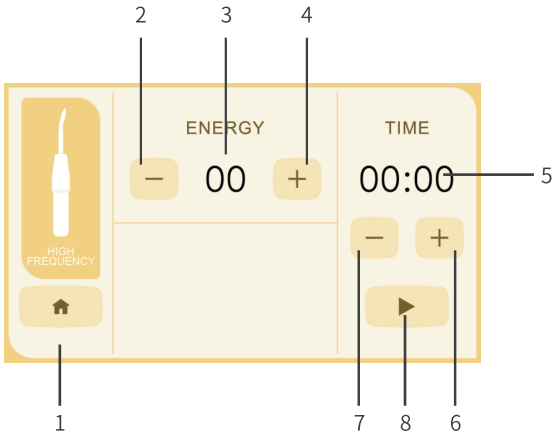

1. Return to main interface
2. Reducing Energy Level
   Energy level range:1~10
3. Energy Level Display
4. Increasing Energy Level
   Energy level range:1~10
5. Working Time Display
   The working time enters the countdown when operating.
   The machine will stop automatically when the time shows zero.
6. Increasing working time
   Working time range:1~30 minutes.
7. Reducing working time
   Working time range:1~30 minutes.
8. Start/Pause

Operation Method

1) Set the output gear and working time, the longest setting can be 30min;
2) Take out the handle and install the electrotherapy tube needed;
3) Press the Start button to enable the function;
4) After the scheduled end, the accessories stop working, take out the electrotherapy tube, clean and wipe the high frequency electrotherapy accessories probe, and then put it back to the accessories shelf.

Note: When holding the high-frequency electrotherapy handle, pay attention to grasp the finger groove, and do not grasp the handle near the probe of the electrotherapy rod attachment to avoid electrical injury.

Recommended using time

① Oily skin: recommended to use 10-15 minutes;
② Dry and dehydrated skin: 3-5 minutes is recommended.

Troubleshooting

1. In case of abnormality, if there is burning smell, stop using immediately and disconnect the power supply.
2. After normal power on, before any function is started, the main unit heats up or shakes obviously for no reason, please disconnect the power immediately and stop using, and contact the maintenance center or supplier to troubleshoot the problem.

The Fault Phenomenon	Reason	Solution
The main machine does not work or crashes.	1. No power.	1. Check whether the fuse in the fuse bin in the power base is burned out. If burn out, replace the corresponding specificat ons of the new fuse. 2. Press the start button, the button will light up after power.
	2. The power cord is not properly inserted or connected.	1. Plug or connect the power cord properly.
Display screen doesn't work.	1. No power.	1. Turn on the power.
The two probes and high-frequency doesn't work	Poor contact of the connectors	Unplug the cable and then reinsert the socket
The output of high-frequency is weak	1. The gear is not set well 2. The electrotherapy stick is aging after long time use	1. Adjust higher gear 2. Replace a new electrotherapy stick
The output of two probes are weak	1. The gear is not set well 2. The temperature is too high, resulting in low output power 3. The probes are damaged	1. Adjust higher gear 2. Use appropriately, do not open for a long time. Do not enable the function when the attachment is not in use. 3. Replace a new probe

If problem is still not solved after above operations, please contact maintenance center or supplier to solve it.



Do not dispose of old or defective appliances in domestic garbage;
this should only be done through public collection points.

▲ Please check the product specification before use. This manual is for reference only. Specifications are subject to change without prior notice.